

Fall Monitors

Whether in hospital or senior living environments, Rondish fall monitors serve as an important element of any comprehensive fall management program.



CordFree NGM-21

1-on-9



-  Monitor up to 9 sensor devices simultaneously
-  Or act as a room controller unit to manage up to 9 beds
-  Sleep mode enables the flexibility to switch off the alarm for low falls risk patients



Palomino NGM-11

1-on-1



-  Unique change pad notification, will alert caregivers prior to expiration of timed sensor pad
-  Works with wetness sensor pad (coming soon) for incontinence care

BedWatcher™

Bed Monitoring System for Fall Prevention

We aim to improve the quality of life for elderly, patients and caregivers with human-centric technology.

Common Features

- Pause mode allows patient to momentarily move off the sensor pad without triggering alarm
- Low battery / signal lost alerts with sensor ID
- Delay settings 0, 2s & 15 mins helps prevent false alarm
- Patented tamper-proof reset option prevent unauthorized cancellation of the alarm
- Nurse call jack allows direct notification to nurse call
- 3 types of tones and 3 levels of volume settings to choose from
- Operates on batteries or AC adaptor




InvisaPad

InvisaPad is designed to work with almost all mattresses, including airflow, foam, spring and extra-thick varieties. It is placed beneath a mattress, so is invisible to a patient. It detects bed presence or absence of a patient without interrupting sleep. Provide patients with the best comfort.



Universal Chair Pad

The pad is 10 by 15 inch with pressure sensors built in. It is placed on or discreetly under a seat cushion to alert caregivers when a patient exits the chair. The pad is designed to be water and incontinence proof.

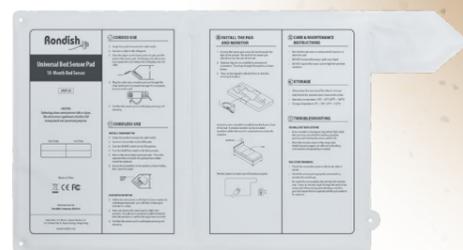

Universal Bed Pad

The pad is 10 by 30 inch with pressure sensors built in. It is placed on or discreetly under the bed sheet to alert caregivers when a patient exits the bed. The pad is designed to be water and incontinence proof.



Early Alert Sensor Pad

Different from the conventional sensor pads, this sensor pad gives an early warning when a patient sits on the edge of the bed, providing a caregiver more time to get to the patient before he/she can make it out of bed on their own. In addition the pad has a delay in the center area to help prevent false alarms. Low false alarms but early alert on leaving the bed for the best solution on the market.





PIR Infrared Sensor

An alternative to sensor pads, this Passive Infrared (PIR) alarm system provides you with simple and dependable motion detection to help manage fall prevention. The PIR alarm has a thin, curtain-like motion detector beam, which only monitors a narrow area beside the bed. When the beam is interrupted, the sensor detects motion and triggers the alarm. As with all contactless devices there are no infection control issue.



CordFree Monitor

The CordFree monitor acts as a room controller unit. Manages up to 9 sensor devices including bed pads, chair pads, floor mats and PIR sensors. Provide you an all-around monitoring. The unit can be set up as a stand-alone solution, or you can configure signals to be received by the central monitoring system.

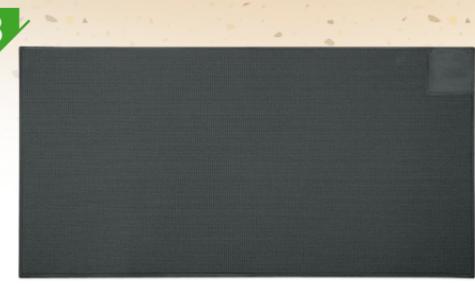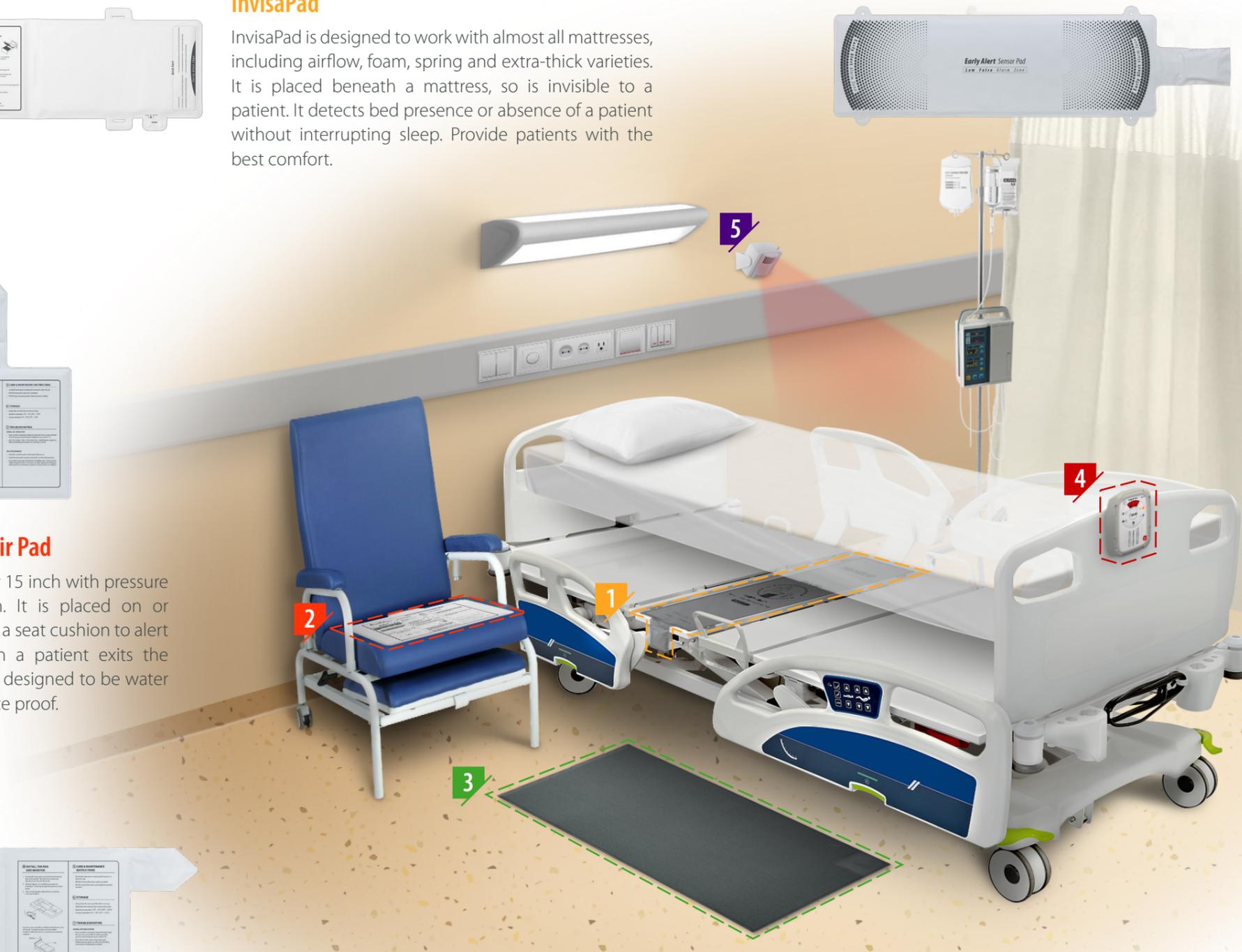

Floor Sensor Mat

Weight-sensing floor mat are placed next to a bed or in a doorway and connected to a monitor. When a patient gets out of the bed or attempts to leave the room and pressure is applied by stepping upon the mat, the monitor will alarm to notify caregivers.

System Components

Cordless Monitor

NGM-21  **CordFree Monitor**
Monitor Up to 9 Sensor Devices • Pause Mode • Sleep Mode • Low Battery Alert • Signal Lost Alert • Tamper-Proof Reset • Delay Setting • 3 Tone Selections • Adjustable Volume • Nurse Call Output • Operates on AC Adaptor or Batteries • Protective Cover • Hanger Strap

NGM-41-FSK  **Nexus Max Monitor**
Monitor Up to 3 Sensor Devices • Pause Mode • Low Battery Alert • Signal Lost Alert • Adjustable Volume • Nurse Call Output • On/Off Switch • Operates on AC Adaptor or Batteries • LoRa Transmission

NGM-41  **Safety First Monitor**
Monitor Up to 3 Sensor Devices • Low Battery Alert • Signal Lost Alert • Adjustable Volume • Nurse Call Output • On/Off Switch • Operates on AC Adaptor or Batteries

Corded Monitor

NGM-11  **Palomino Monitor**
Change Pad Notification • Incontinence Alert • Pause Mode • Low Battery Alert • Input Lost Alert • Tamper-Proof Reset • Delay Setting • 3 Tone Selections • Adjustable Volume • Nurse Call Output • Operates on AC Adaptor or Batteries • Protective Cover • Hanger Strap

BMS-03  **Secure Reset Safemon**
Low Battery Alert • Input Lost Alert • Tamper-Proof Reset • Tone Selection • Adjustable Volume • Nurse Call Output • Nurse Call Plug Removal Alert • Operates on AC Adaptor or Batteries • Protective Cover • Hanger Strap

Sensor Devices

UFBP-20  **Universal Bed Pad**
Corded & Cordless 2-in-1 • Foldable Design • 30" x 20" (76 x 51cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)

WUMP-40  **InvisaPad**
Works Underneath the Mattress • Auto Calibration Technology • Foldable Design • 32" x 9.5" (81 x 24cm) • Optional Early-Alert feature • Built-in Transmitter

UMBP-10-U  **Universal Under-Mattress Bed Pad**
Works Beneath Most Mattresses • Corded & Cordless 2-in-1 • 33" x 10" (84 x 25cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)

WPAS-10  **Early-Alert Sensor Bed Pad**
Sends Alarm BEFORE Patient Has Left The Bed • Multi-Zone Detection Design • 30"x10" (76 x 25 cm) • Water Resistant • Built-in Transmitter

UCP-10  **Universal Chair Pad**
Corded & Cordless 2-in-1 • 15" x 10" (38 x 25cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)

TM-12U  **Transmitter**
Delay Setting • Switch to Change Between Bed & Chair Type Alarms • Built-in Batteries • Low Battery Indication • Mini USB Plug

UFFMT-36  **Universal One-fold Floor Mat**
Corded & Cordless 2-in-1 • Foldable Design • 36" x 24" (90 x 60cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)

TM-12BU  **Transmitter**
Delay Setting • Switch to Change Between Bed & Chair Type Alarms • Replaceable Batteries • Low Battery Indication • Mini USB Plug

UFMT-10  **Universal Heavy-duty Floor Mat**
Corded & Cordless 2-in-1 • 48" x 24" (120 x 60cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)

USBC-0.8M  **Breakaway Cable**
0.8m Mini USB Breakaway Cable with RJ10 Plug

PIR-21  **Passive Infrared Sensor**
LED Status Indicator • Internal Alarm Or Signal Transmit Options • Pause Mode • Easy Reset Button • Adjustable Detection Range • 3xAA Battery Operated • Low Battery Alert • DC (6V) Input Option

USBC-2M  **Breakaway Cable**
2m Mini USB Breakaway Cable with RJ10 Plug

TXP-11(V2)  **Wristband Transmitter**
Man-down Detection • Active Emergency Call Button • Waterproof to IP65 • Replaceable Battery • Comes With Faux Leather Watch Strap



We aim to improve the quality of life for elderly, patients and caregivers with human-centric technology.



(852) 2541 6685

(852) 2541 7411

sales@rondish.com

www.rondish.com

facebook/RondishCare

youtube.com/@RondishCare