



■ Nexus Master Station

NEXUS

Nexus is a **Wireless Nurse Call System** that is designed based on UL 2560.

The solution is developed from ground up to provide a **Reliable & Price Competitive** wireless connection. By employing the latest **LoRa (Long Range) Technology**, Nexus creates **Larger Coverage** for broader and more extensive communication. The system now includes **Fall Management & Wander Management Solutions** to provide greater performance capabilities.



Touchscreen Master Station

- Color Coded Touchscreen Display
- Calls Ordered By Priority
- Audio Call Levels Indication
- Low Battery And Device Lost Alerts
- Event Log and Response Time Tracking
- Password Protected
- Options for Wall-mount / Desk-mount

NEXUS Key Features & Benefits :

- 
Simple & Quick Self-Installation
 Console, call points and sensor devices need no wiring. Eliminates installation complexity
- 
Easy-To-Use
 Intuitive UI design, straightforward to user and no specialized training required
- 
Extensive Coverage
 Open air range in excess of one mile
- 
Scalable & Flexible
 Unlimited capacity of call points and wide range of sensor devices that can cover an extensive physical area
- 
Easy Deployment
 Fully functional nurse call system can be deployed within a matter of minutes
- 
Fall / Wander Management Integration
 Integrates seamlessly with the fall / wander management solutions. Direct alerts to the nurse call system for central monitoring and response. Local alerts can be disabled for silent alarming
- 
OverDoor Indicator
 Visual and audible alarm, when acts as a signal repeater for extending coverage
- 
Fault Supervision
 Alerts for product faults, out of range and low battery conditions
- 
Event Log
 Every event that occurs on the system is collected. Data can be exported for analysis and reporting
- 
Reliable Backups & Rapid Recovery Capabilities
 Instant recovery solution in the event of an unexpected hardware or system failure
- 
Messaging Solutions Available
 Instant notification sent to staff members via in-house pagers or smartphones
- 
Real-time Multiple Display Panels
 Staff do not have to travel far to see the location of the call

Assisted Living Facility



Fall Management

Bedhead Call Point

Nurse Station

Wearable Device


Wander Management

Bathroom / WC Safety



Fall Monitor

The CordFree monitor acts as a room controller unit. Manages up to 9 sensors and supports a wide range of sensor devices

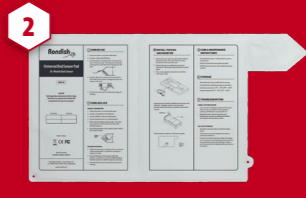


Sensor Pad

Various types & sizes of sensor pads, include above or under-mattress bed pads, chair pads, floor mats are available




Portable Handset

A portable nurse call handset can be fastened on a bed rail. It provides the most affordable solution for bedding needs

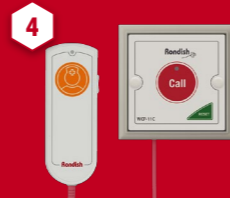

Bedhead Call Point

Call & Reset 2-in-1 bedhead unit is completely wireless. Battery operated with long battery life expectancy. 4P4C jack for a nurse call handset

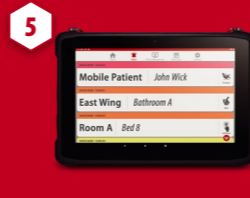

Master Station

A touchscreen Master Station indicates calls in order of priority. Can either stand on a desk or wall-mounted



Wristband Transmitter

The Wristband Transmitter with man-down detection and emergency call button allow patients to call for assistance. Pendant is available for option.



Door Strip Sensor & Wristband Transmitter

Provides audible and visual alert when patient wearing wristband transmitter move within the range of protected doors



OverDoor Indicator

Provides visual and audible alarm, while acting as a signal repeater for extending coverage



Bathroom Call Point






Wireless Call & Reset 2-in-1 button is waterproof to IP65. It comes with a pull cord to ensure accessibility from the toilet or floor

Multiple Options • Proven Performance


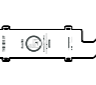





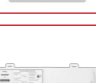

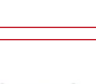
System Components

 <p>NDP-21(V2)</p>	<p>Master Station 10.1" Rugged Tablet PC • With Nexus Software Installed • IP65 • 8000mAh Battery</p>	 <p>ILB-21-FSK</p>	<p>Overdoor Indicator Dual Colored Indicator Light • Buzzer With Adjustable Volume • Signal Repeater • Battery Back-up In Case Of Power Outage</p>
--	--	--	---



Nurse Call Components

 <p>WCP-11C</p>	<p>Bedhead Call Point Call & Reset Buttons 2-in-1 • Call Reassurance LED • 4P4C Connection For Nurse Call Handset NEC-12(M) • Replaceable Batteries • Long Battery Life</p>	 <p>NEC-12(M)</p>	<p>Nurse Call Handset Easy To Press Big Amber Call Button • Call Reassurance LED • Standard Cable With USB Break Point & 4P4C Connector • Comes With Wall-mounted Bracket • Compatible With WCP-11C</p>
 <p>WCP-11</p>	<p>Bathroom / WC Call Point Call & Reset Buttons 2-in-1 • Comes With Removable Pull Cord • IP65 Waterproof Rating • Replaceable Batteries • Long Battery Life</p>	 <p>NEC-05(M)</p>	<p>Portable Handset Easy To Press Big Amber Call Button • Grey Button To Reset • Call Reassurance LED • 2xAAA Batteries Operated • Comes With Clip & Plastic Wall-mounted Bracket for Storage</p>
 <p>TXP-11(V2)</p>	<p>Wristband Transmitter Man-down Detection • Active Emergency Call Button • Waterproof to IP65 • Replaceable Battery • Comes With Faux Leather Watch Strap</p>		

Fall Management Components

 <p>NGM-21</p>	<p>CordFree Monitor Monitor Up to 9 Sensor Devices • Pause Mode • Sleep Mode • Low Battery Alert • Signal Lost Alert • Tamper-Proof Reset • Delay Setting • 3 Tone Selections • Adjustable Volume • Nurse Call Output • Operates On AC Adaptor or Batteries • Protective Cover • Hanger Strap</p>	 <p>UMBP-10-U</p>	<p>Universal Under-Mattress Bed Pad Works Beneath Most Mattresses • Corded & Cordless 2-in-1 • 33" x 10" (84 x 25cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)</p>
 <p>NGM-41-FSK</p>	<p>Nexus Max Monitor Monitor Up to 3 Sensor Devices • Pause Mode • Low Battery Alert • Signal Lost Alert • Adjustable Volume • Nurse Call Output • Operates On AC Adaptor or Batteries</p>	 <p>UFBP-20</p>	<p>Universal Bed Pad Corded & Cordless 2-in-1 • Foldable Design • 30" x 20" (76 x 51cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)</p>
 <p>UFMT-10</p>	<p>Universal Heavy-duty Floor Mat Corded & Cordless 2-in-1 • 48" x 24" (120 x 60cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)</p>	 <p>UCP-10</p>	<p>Universal Chair Pad Corded & Cordless 2-in-1 • 15" x 10" (38 x 25cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)</p>
 <p>UFFMT-36</p>	<p>Universal One-fold Floor Mat Corded & Cordless 2-in-1 • Foldable Design • 36" x 24" (90 x 60cm) • Water Resistant • Mini USB Plug for Breakaway Cable or Rondish Transmitter (TM-12U / TM-12BU)</p>	 <p>WUMP-40</p>	<p>InvisaPad Works Underneath the Mattress • Auto Calibration Technology • Foldable Design • 32" x 9.5" (81 x 24cm) • Optional Early-Alert feature • Build-in Transmitter</p>
 <p>PIR-21</p>	<p>Passive Infrared Sensor LED Status Indicator • Internal Alarm Or Signal Transmit Options • Pause Mode • Easy Reset Button • Adjustable Detection Range • 3xAA Battery Operated • Low Battery Alert • DC (6V) Input Option</p>	 <p>WPAS-10</p>	<p>Early-Alert Sensor Bed Pad Sends Alarm BEFORE Patient Has Left The Bed • Multi-Zone Detection Design • 30"x10" (76 x 25 cm) • Water Resistant • Built-in Transmitter</p>

Wander Management Components

 <p>DMS-21</p>	<p>Door Strip Sensor Visual & Audible Alert • Adjustable Sensing Range • Automatic Door Lock Option • Customize ID • Low Battery Alerts • Fault-checked</p>	 <p>TXP-21</p>	<p>Wristband Transmitter Unique tag ID • Waterproof to IP67 • Low Battery Alert • Comes With Faux Leather Watch Strap</p>
--	--	--	--



We aim to improve the quality of life for elderly, patients and caregivers with human-centric technology.

(852) 2541 6685

sales@rondish.com

facebook/RondishCare

(852) 2541 7411

www.rondish.com

youtube.com/@RondishCare