

HINGMED

DBP-01P/DBP-01HP

Clinical Automatic BP Monitor




Make everyone's blood pressure healthier


Advantages & Features



- Pulse wave anti-interference technology, accurate measurement.
- Pulse wave intelligent booster measurement, fast & comfortable.
- Ultra short measurement time, 20 seconds.

- Wide range arm circumference, automatic suit different arm.
- The sleeve can be adjusted up and down, left and right.
- Infrared detection for pre-readiness.

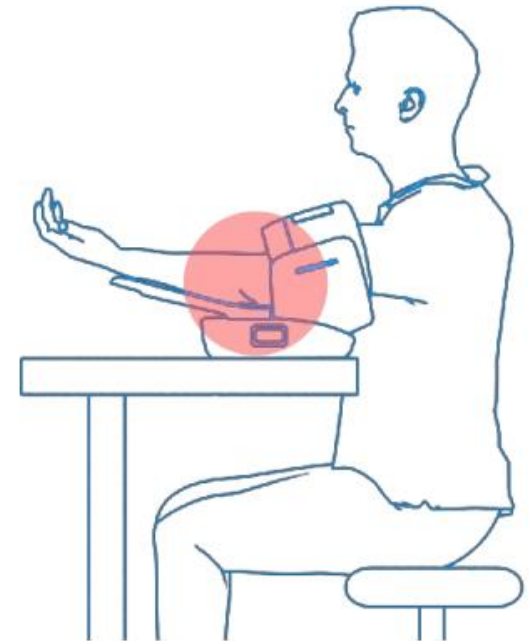

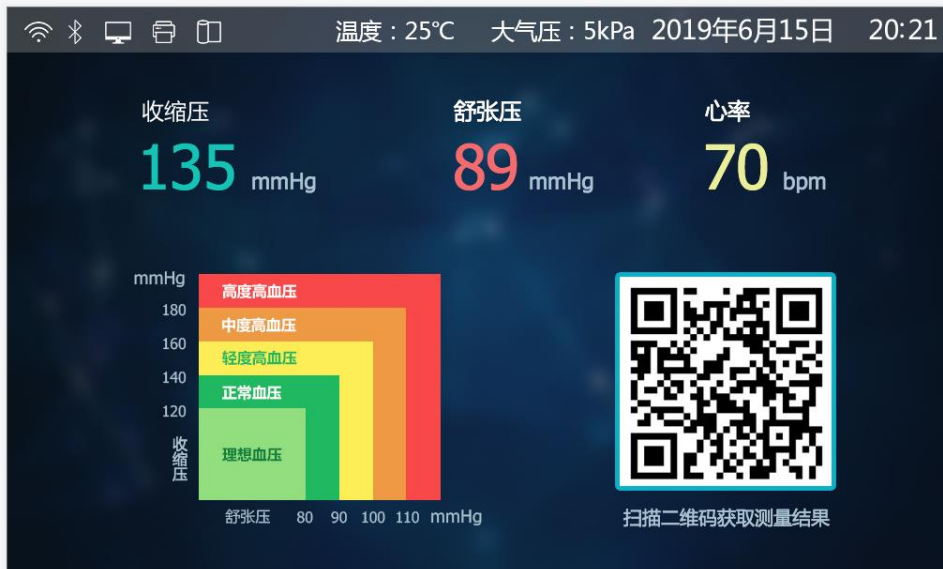

- 8" colorful display, support video playing. (DBP-01H)
- Android system support WiFi. (DPB-01H)
- Pressure & Temperature indication


- Various connection port:
USB, RS232, ETHERNET, SD card.





- Modular Design.
- Easy maintenance.

- Clinical automatic BP module.
- Small size, easy connection.

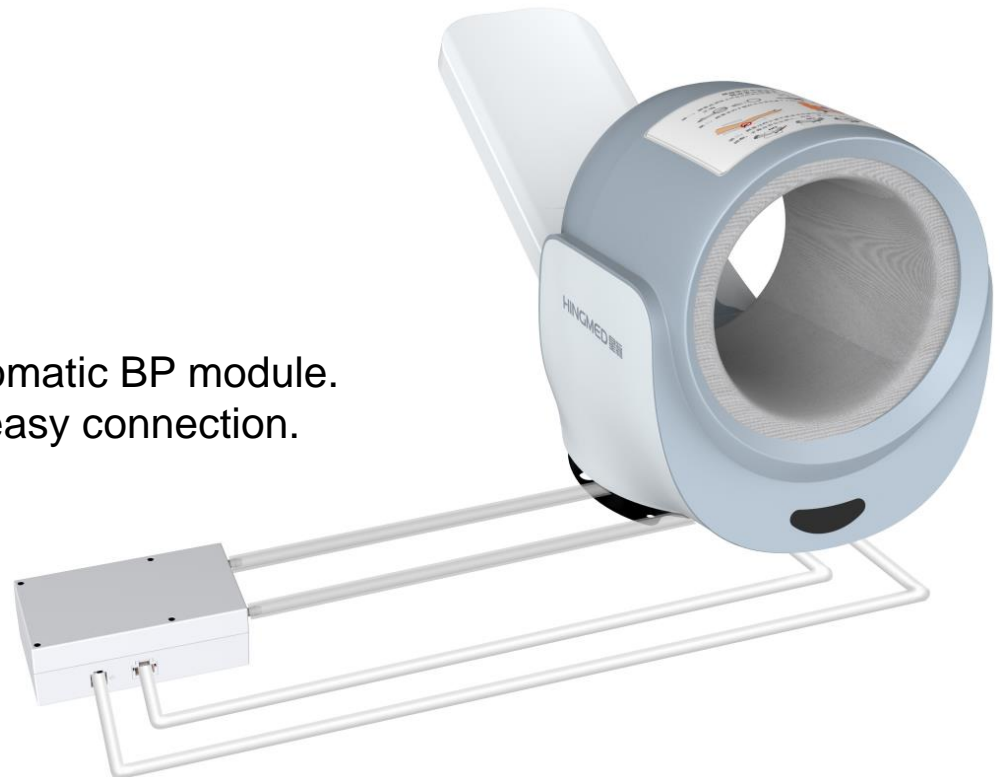

Specification

Name	Clinical Automatic BP Monitor
Model	DBP-01/DBP-01P/DBP-01H/DBP-01HP
Technology	Intelligent Inflation Oscillometric
Systolic	40-260 mmHg
Diastolic	20-210 mmHg
Pulse Rate	40-200 bpm
Accuracy	±3 mmHg
Calibration	At least once per two years
Safety	Valve open automatically for deflation upon max inflation pressure 290mmHg.Pushing “Stop”button upon emergency,deflating quickly.
Memory	9999 readings
Power	110-240V AC,50/60hZ
Printer	thermal printing,Paper size58mm
Cuff Circumference	17-42cm
Voice	Voice introduction of operation & reading,Language Option:Chinese,English.On/Off Option

Size	310mm*476mm*300mm
Weight	6.5kg
Electric Shock Protection	Class I
Working Condition	Tem:+5°C--+40°C, RH: 10%-95%, Air Pressure 70-106Kpa.
Store Condition	Tem: -2°C--+55°C, RH: ≤95%, Air Pressure 70-106Kpa

Company: Shenzhen Hingmed Instruments Co., Ltd.

Address: 4th Floor, Zhonghangfeixiang Building, NO. 371, Guangshen Road, Baoan District, Shenzhen, People's Republic of China

Tel: +86 755 23003372 Post code: 518102